

FINAL REPORT

TIME KILL STUDY

PROCEDURE NO. STP0158 REV 01  
PROTOCOL DETAIL SHEET NO. 200901442 REV 01

LABORATORY NO. 474527E

PREPARED FOR:

KEITH MOELLER  
AMERICAN BIOTECH LABS  
80 W. CANYON CREST RD.  
ALPINE UT 84004

SUBMITTED BY:

NELSON LABORATORIES, INC.  
6280 S. REDWOOD RD.  
SALT LAKE CITY UT 84123-6600  
801-290-7500



NELSON LABORATORIES, INC.

QAU AUDIT STATEMENT

☒ USFDA (21 CFR PART 58)

☐ USEPA (40 CFR PART 160)

TIME KILL STUDY

LABORATORY NO. 474527E

1. The test was conducted in accordance with the USFDA or USEPA Regulations as noted above.
2. In accordance with the Good Laboratory Practice Regulations, the Culture Prep phase(s) of this study was inspected by the Quality Assurance Unit on: 13 May 2009. The findings of the inspection(s) were reported to the Study Director and to Management on: 09 Jun 2009.
3. The Quality Assurance Unit has reviewed this report and has determined that the methods and standard testing procedures are accurately described, and that the reported results accurately reflect the raw data.
4. The name of the study director, the names of other scientists or professionals, and the names of all supervisory personnel, involved in the study:

Mike Neilson  
Peter Croci

Dr. Jerry Nelson  
Jeff Hills

QUALITY ASSURANCE:

Brent M. Smith

DATE: 16 Jun 2009



## TIME KILL STUDY

LABORATORY NUMBER:	474527E
PROCEDURE NUMBER:	STP0158 REV 01
PROTOCOL DETAIL SHEET NUMBER:	200901442 REV 01
SAMPLE SOURCE:	American Biotech Labs
SAMPLE IDENTIFICATION:	Refer to Tables 1-2
DEVIATIONS:	None
PROTOCOL APPROVAL DATE:	08 May 2009
SAMPLE RECEIVED DATE:	07 May 2009
LAB PHASE START DATE:	08 May 2009
LAB PHASE COMPLETION DATE:	09 Jun 2009
REPORT ISSUE DATE:	11 Jun 2009

### INTRODUCTION:

This report describes the procedure for the evaluation of products for anti-microbial activity against selected organisms at representative contact times. Products are evaluated in a liquid matrix. The test organisms and contact times are chosen by the sponsor. This is a quantitative test that allows the determination of the amount of organism reduction at pre-determined intervals.

### ACCEPTANCE CRITERIA:

Negative controls must be negative for growth. Positive controls must be positive for growth. Neutralization must be confirmed at  $\geq 70\%$ . Specific criteria for pass/fail of the sample must be determined by the Sponsor.

### INOCULUM PREPARATION:

Plates of brain heart infusion agar (BHIA) media were inoculated with stock cultures of bacteria and incubated at 30-35°C for 18-24 hours. Growth was harvested from the surface using a sterile bent glass rod and purified water solution (PURW).

Where necessary, culture suspensions were adjusted for the test procedure with purified water solution (PURW) to approximately  $(\sim) 10^8$  colony forming units (CFU)/mL using visual turbidity.



American Biotech Labs  
Lab Number 474527E

Time Kill Study

#### SAMPLE PREPARATION:

Samples were prepared according to the product label or sponsor instructions and were tested without any additional manipulation or dilution.

#### NEUTRALIZATION:

A 1 mL aliquot of the test sample was mixed with 9 mL of Dey-Engley neutralizer broth (DEYE). An additional tube of 10 mL of DEYE was prepared as a titer control. The tubes were inoculated with 0.1 mL of a test organism suspension diluted to approximately  $\leq 10,000$  CFU/mL, then mixed thoroughly. Aliquots from each tube were plated in triplicate onto BHIA and incubated as described in the test procedure.

The neutralizer is considered effective for the product if the number of CFU on the test plates is  $\geq 70\%$  of the control plates.

#### CONTROLS:

Positive control tubes containing 10 mL physiological saline (PHSS) were prepared. 0.1 mL of the test organism was added to each tube. The positive control was assayed at 0 and 24 hour.

The negative control consisted of plating sterile aliquots of applicable liquid media in triplicate and incubating as described in the test procedure.


#### TEST PROCEDURE:

Tubes containing 10 mL of each test sample were prepared and inoculated with 0.1 mL of the test organism to yield  $\sim 10^8$  CFU/mL. The samples were mixed thoroughly.

At 1 hour and 24 hours of exposure 1.0 mL aliquots of the test suspension were removed and added to 9 mL of neutralizer. The tubes were mixed thoroughly. Serial dilutions were made in the appropriate neutralizer and assayed using a standard spread plate method.

All plating was performed in triplicate. Bacterial samples were plated onto BHIA and incubated at 30-35°C for 2-5 days.



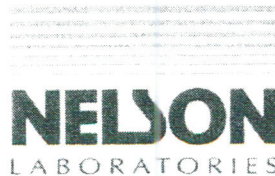

American Biotech Labs  
Lab Number 474527E

Time Kill Study

#### CALCULATIONS:

The log reduction values were calculated using the following formula:

$$\log \text{reduction} = \log U - \log C$$

Where U = Average positive control titer

Where C = Average recovered counts

The percent reduction values were calculated using the following formula:

$$\% \text{ reduction} = 1 - \frac{1}{10^{(\log \text{ reduction})}} \times 100\%$$

The percent neutralization is obtained according to the following equation:

$$\% \text{ Neutralization} = \frac{\text{Average Sample Counts/Plate}}{\text{Average Control Counts/Plate}}$$

#### RESULTS:


All testing was performed in accordance with Nelson Laboratories, Inc. (NLI) protocol detail sheet #200901442 REV 01.

The average sample and control counts, percent reductions and log<sub>10</sub> reductions for the samples can be found in Tables 1. Values are considered approximate (~) when plate counts were outside of the statistically accurate range of 25-250 CFU/plate for bacteria and yeast and 8-80 CFU/plate for mold. Less than symbols (<) are applied to recovery values where no CFU were observed on the plates. This denotes the limit of detection for the test. Neutralization results are summarized in Table 2 with all samples demonstrating ≥ 70% recovery.

Testing met the acceptance criteria previously described in this report.

#### CONCLUSION:

Interpretation of the data is the responsibility of the sponsor and no conclusion can be made by Nelson Laboratories, Inc. (NLI).

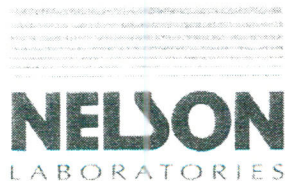

American Biotech Labs  
Lab Number 474527E

Time Kill Study

DATA DISPOSITION:

The raw data and final report from this study are archived at NLI or an approved off-site location.

STATEMENT OF UNCERTAINTY

If applicable, a statement of uncertainty is available to sponsors upon request.

A handwritten signature in dark ink, appearing to read "Peter Croci", written over a horizontal line.

Peter Croci, B.A.  
Study Director

A handwritten date "11 Jun 2009" in dark ink, written over a horizontal line.

Study Completion Date

PC/es

**NELSON**  
LABORATORIES

American Biotech Labs  
Lab Number 474527E

Time Kill Study

TABLE 1. Results  
Sample Identification: ASAP Wound Dressing Gel 050509-1

ORGANISM	EXPOSURE INTERVALS	AVERAGE CONTROL TITER (CFU/mL)	PERCENT REDUCTION (%)	LOG <sub>10</sub> REDUCTION
Vancomycin-resistant <i>Enterococcus faecalis</i> , ATCC #51575	1 Hour	4.8 x 10 <sup>6</sup>	99.56	2.35
	24 Hour		>99.99958	>5.38
	24 Hour Control		5	0.02

TABLE 2. Neutralization

SAMPLE IDENTIFICATION	ORGANISM	AVERAGE CONTROL COUNTS (CFU)	AVERAGE SAMPLE COUNTS (CFU)	PERCENT NEUTRALIZATION (%)
ASAP Wound Dressing Gel 050509-1	Vancomycin-resistant <i>E. faecalis</i>	53	45	85



American Biotech Labs  
Lab Number 474527E

Time Kill Study

All reports and letters issued by Nelson Laboratories, Inc. are for the exclusive use of the sponsor to whom they are addressed. These results relate only to the samples tested. Reports may not be reproduced except in their entirety. No quotations from reports or use of the corporate name is permitted except as expressly authorized by Nelson Laboratories, Inc. in writing. The significance of any data is subject to the adequacy and representative character of the samples tendered for testing. Nelson Laboratories, Inc. warrants that all tests are performed in accordance with established laboratory procedures and standards. Nelson Laboratories, Inc. makes no other warranties of any kind, express or implied. Nelson Laboratories, Inc. expressly states that it makes no representation or warranty regarding the adequacy of the samples tendered for testing for any specific use of application, that determination being the sole responsibility of the sponsor. Nelson Laboratories' liability for any loss or damage resulting from its actions or failure to act shall not exceed the cost of tests performed, and it shall not be liable for any incidental or consequential damages. Nelson Laboratories, Inc. has a ten year record retention policy.